

INSPIRATION

Radiant Heat Warmers

NEXT ONLY TO MOTHER'S CARE

With early warning of
Impending Cold Stress Onset (Optional)

Key Features

- On Board Temperature Calibration Facility with Feather Touch Keys
- Entire Warmer Height Up-down Facility with Motorized Remote Control (Optional)
- Baby Tray Up-down (Trendelenburg) Facility (Optional with Motorized Remote Control)
- Swiveling 3 Drawers



Micro Computer Temperature Controller

- PID (proportionate integrated derivative) based micro processor technology
- Dual monitoring of Central and Peripheral temp. with early warning of Cold Stress onset
- User friendly on-board temperature calibration with feather touch keys
- User friendly feather touch control
- Key-pad lock facility
- Safety alarms along with high temp. cutoff
- Watch Dog facility

Heater Box


- Ceramic Infra Red Heater
- Swivelling source box to take X-ray
- Stainless steel parabolic reflector with heater grill
- LED Examination Lamp

Baby Bassinet

- Easy access of an Infant from all sides due to clear acrylic collapsible side panels
- Provision for positioning X-Ray cassette
- Thick cushion mattress

Essential Attachments

- Instrument Tray
- Swiveling 3 Drawers
- Pole for Infusion Pump & IV Stand

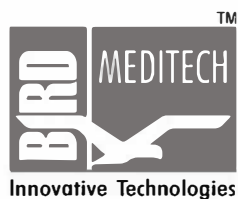

Bird Meditech
ISO 13485:2012 CERTIFIED COMPANY

www.birdmeditech.com



OPTIONAL ATTCHMENT



- ❖ Masimo SPO₂
- ❖ Inline Weighing
- ❖ T-Piece Resuscitator
- ❖ Neonatal Suction
- ❖ Heater Source remains Focused on Baby even in Swivel condition

Displays

- Skin temperature
- Air temperature
- Peripheral temperature (Optional)
- Heater output (%)
- Set temperature
- Control mode
- APGAR Timer

Manual Mode

- Time programmable
- Heater output adjustable
- Over temp. cutoff.
- Timer / APGAR Timer
- Baby attending alarm
- Mute time programmable
- Audio alarm with mute facility

Power Consumption : 800 Watt

Auto Mode

Safety alarm messages on LCD

Skin Mode

"Skin Temp. High" (>0.5°C of set temp.)
"Skin Temp. Low" (<0.5°C of set temp.)
"Skin Temp. Over Range" (>38°C)
"Skin Temp. Over Range" (<34°C)
"Skin Probe Fail"

Cold Stress Alarm (Optional)

Peripheral Temp. 1.5°C < Skin Temp.

Air Mode

"Air Temp. High" (>1.5°C of set temp.)
"Air Temp. Low" (<3.0°C of set temp.)
"Air Temp. Over Range" (>39°C)
"Air Probe Fail"

Safety high temp. cutoff at

Skin temp. >38°C
Air temp. >39°C
Power Fail
Heater Fail (Optional)

Available at:

Xavier Med Pvt. Ltd.

1-10-51/2, Road 6B, Dwarakapuram,
Dilsukhnagar, Hyderabad - 500060

www.buyxavier.com

Customer Care: 8008895680,
Service: 8008895079